

*cryo*SPHERE®

- Enhanced Flexibility
- Uniform Application
- Engineered for Intercostal Nerve Cryoablation

Indications

AtriCure's cryoICE® cryoSPHERE® cryoablation probes are sterile, single use devices intended for use in blocking pain by temporarily ablating intercostal peripheral nerves.

EU Indications: AtriCure's cryoICE cryoSPHERE cryoablation probes are sterile, single use devices intended for use in blocking pain by temporarily ablating intercostal peripheral nerves.

Please review the Instructions for Use for a complete listing of contraindications, warnings, precautions and potential adverse events prior to using these devices.

ATRICURE, INC.

7555 Innovation Way
Mason, Ohio 45040 USA
Customer Service: +1 (866) 349-2342
CustomerService@AtriCure.com
www.AtriCure.com

ATRICURE EUROPE B.V.

De Entree 260
1101 EE Amsterdam
Netherlands
Tel: +31 (0) 20-7005560
Fax: +31 (0) 20-7005561
SalesSupportEU@AtriCure.com
Europe.AtriCure.com

PM-EU-1697D-0227-G


*cryo*SPHERE®
A non-pharmacological/
non-opioid solution that
can provide effective post-
operative pain management
for several months



AtriCure

AtriCure


A non-pharmacological/non-opioid solution that can provide effective post-operative pain management for several months

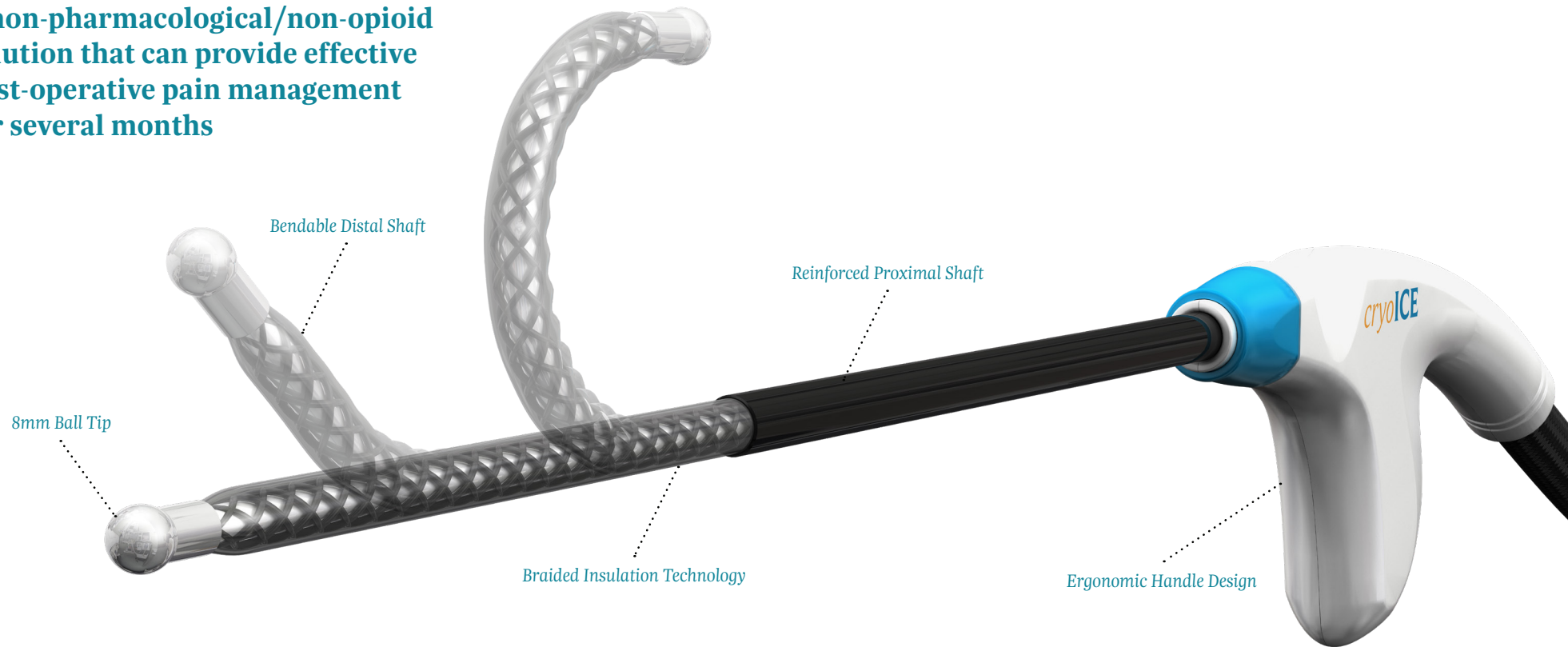

Features

8mm Ball Tip	Designed for minimally invasive applications and placement on intercostal nerve
Bendable Distal Shaft	Bendable distal shaft up to 180 degrees
Braided Insulation	Technology that is designed to maintain a temperature warmer than -10°C in tissue
Reinforced Proximal Shaft	Provides rigidity at the level of the incision points
Ergonomic Handle Design	Provides ease of use and comfort during cryoablation
Active Defrost	Allows probe to be easily removed from cryoablation site

Malleable Cryoablation Probe

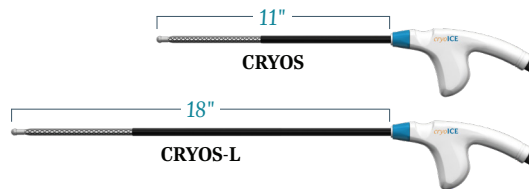

cryoICE BOX®- V6

ACM1 (US)



cryoICE Cart & Upgrade Kits

ASCS (cart with double tank holder, US only)

ASCC (single tank holder upgrade)

ASCC2 (double tank holder upgrade)

20 lbs N₂O tanks not provided by AtriCure



Foot Switch (optional)

For use with ACM1/ACM2
FTSWITCH

