


Isolator® Synergy™ Access®

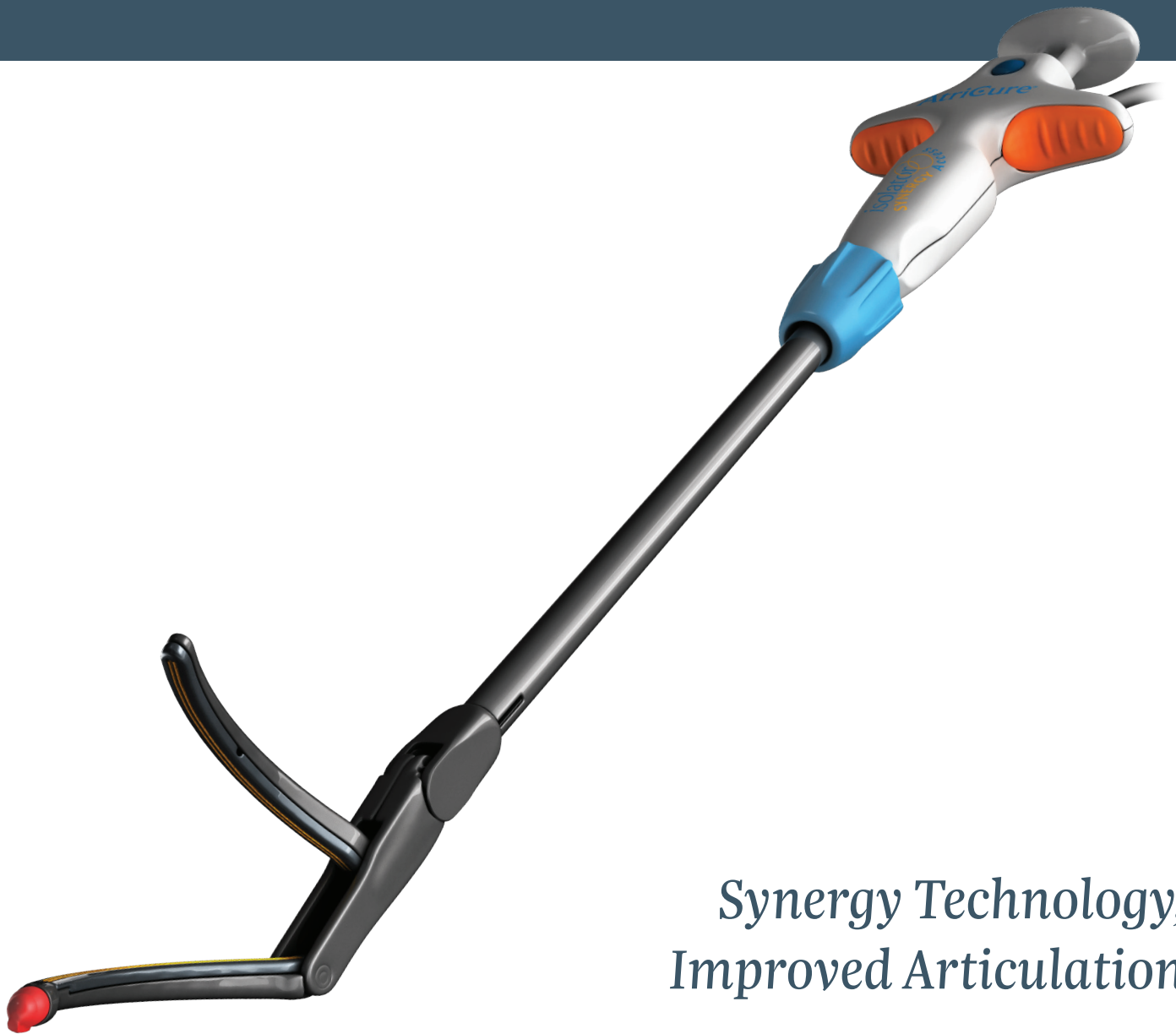


Synergy Technology,
Improved Articulation

AtriCure

Isolator Synergy Access

Isolator Synergy Science

Ground-breaking Synergy Science uses two pair of alternating bipolar electrodes to maintain resistive heating between the four poles. An initial lesion is formed where the RF energy fields overlap and continues until a robust, column shaped lesion is produced.

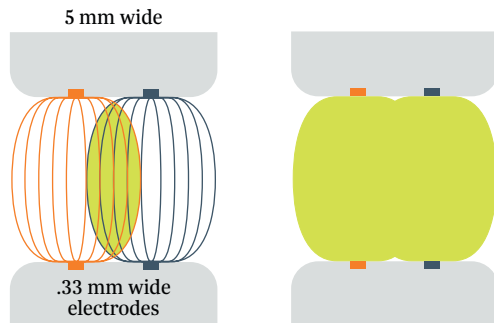

Consistent Tissue Compression

Minimal Tip Splay (<0.06")

+ 0% Jaw Deflection

= True Parallel Closure




ATRICURE EUROPE B.V.
De Entree 260
1101 EE Amsterdam
the Netherlands
Tel.: +31 (0) 20-7005560
Fax: +31 (0) 20-7005561
SalesSupportEU@AtriCure.com
Europe.AtriCure.com

PM-EU-2365B-1226-G

Jaw Articulation

The Synergy Access jaw rotates 30 degrees upward and 30 degrees downward to enable ease of access in even the most difficult anatomy.

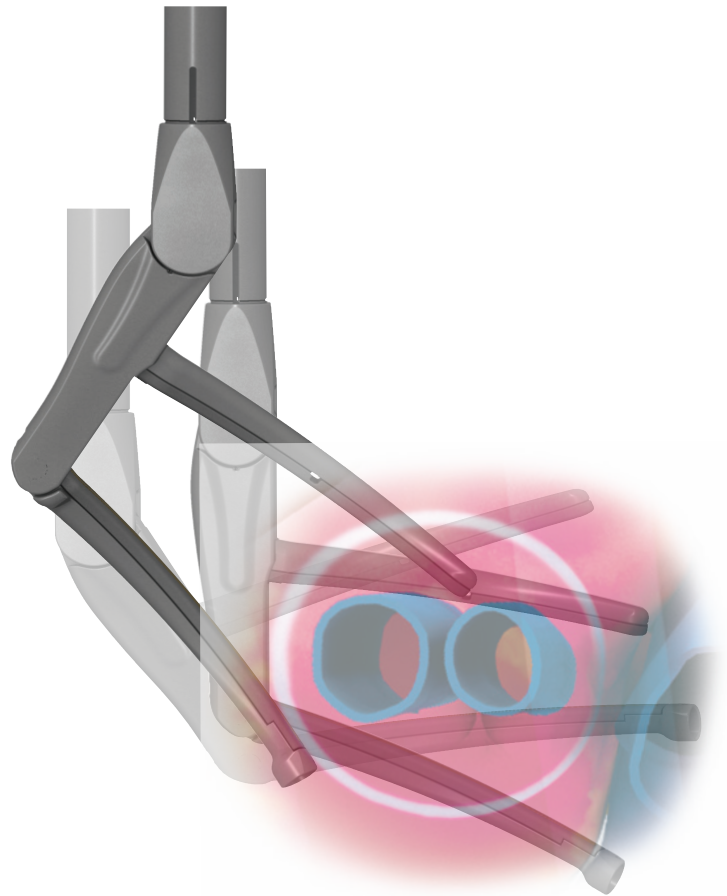

Isolator Synergy Access

Product Code

EMT1

*Includes Glidepath Tape

The ATRICURE Bipolar (Transpolar) System is indicated for ablation and coagulation of soft tissue in general, ENT, thoracic, urological, gynecological surgical procedures and ablation of cardiac tissue during surgery including pulmonary vein isolation and atrial connecting lesions for the Maze procedure for the treatment of cardiac arrhythmias, including atrial fibrillation.

Please review the Instructions for Use for a complete listing of contraindications, warnings, precautions and potential adverse events prior to using these devices.

AtriCure