

AtriClip® Devices

Reduce the Risk of Stroke*

>98%
successful LAA
exclusion¹

AtriCure's AtriClip products are the most widely sold Left Atrial Appendage (LAA) management devices worldwide.

Epicardial Exclusion

- Implant is not in the blood stream
- Ischemic injury electrically isolates the LAA
- Atrophy and resorption of the LAA

Dynamic Closing Force

Continuous closing force maintains LAA exclusion throughout changes to the tissue caused by ischemia

Parallel / Linear Closing


Minimized occurrence of tissue folds with optimal apposition of tissue along long axis of LAA ostia

Tissue Compression / Atraumatic

No cutting, non-piercing, and atraumatic compression reduces risk of tissue tearing and bleeding

Stroke* Risk Reduction

Indicated for use in patients at high risk of thromboembolism for whom left atrial appendage exclusion is warranted

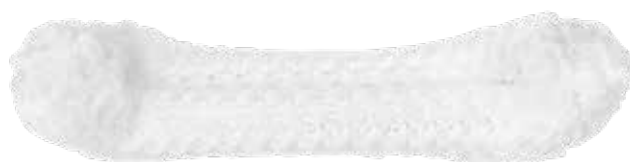

Class I Recommendation
ESC/EACTS 2024 Guidelines

*Originated from the LAA

The Closed-Ended AtriClip Devices

The Open-Ended AtriClip Devices

Open Access

AtriClip FLEX® device

Features

- Plunger grip
- 6 cm flexible shaft
- Available in four sizes



Minimal Access

AtriClip PRO® device

Features

- Quick deploy feature
- 25cm rigide shaft
- ±30° omnidirectional head articulation with lock
- Includes selection guide (CGG)



AtriClip PRO2® device

Features

- Hoopless end effector
- Active articulation levers
- Small diameter end effector (12 mm)



Open Access

AtriClip FLEX•V® device

Features

- Suture-less AtriClip Deployment
- Reduced fatigue Clip Opening Lever
- One handed application



AtriClip FLEX-Mini® device

Features

- 60% smaller profile to minimize interaction with adjacent structures and enhance LAA visualization
- One-handed deployment
- Reduced-fatigue clip opening lever
- Sutureless quick release




Minimal Access

AtriClip PRO•V® device

Features

- 12 mm port compatibility
- Tip-first closure
- Opened-ended design clip

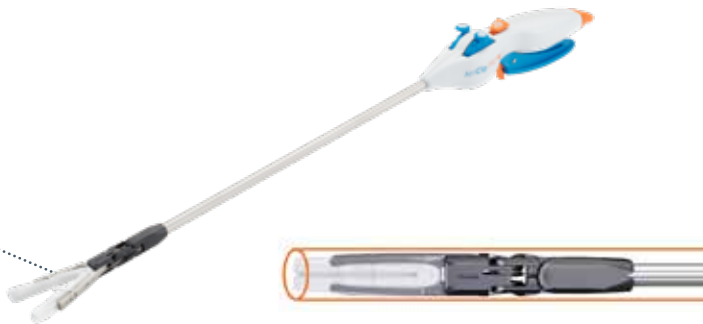

AtriClip PRO-Mini® device

Features

- 60% smaller profile to minimize interaction with adjacent structures and enhance LAA visualization
- 12 mm port compatibility
- Hoopless end effector
- Active articulation for precise positioning



AtriClip FLEX

Device	Size
ACH235	35 mm
ACH240	40 mm
ACH245	45 mm
ACH250	50 mm

AtriClip FLEX•V

Device	Size
ACHV35	35 mm
ACHV40	40 mm
ACHV45	45 mm
ACHV50	50 mm

AtriClip PRO•V

Device	Size
PROV35	35 mm
PROV40	40 mm
PROV45	45 mm
PROV50	50 mm

AtriClip PRO

Device	Size
PRO135	35 mm
PRO140	40 mm
PRO145	45 mm
PRO150	50 mm

AtriClip FLEX-Mini

Device	Size
ACHM35	35 mm
ACHM40	40 mm
ACHM45	45 mm
ACHM50	50 mm

AtriClip PRO-Mini

Device	Size
PROM35	35 mm
PROM40	40 mm
PROM45	45 mm
PROM50	50 mm

AtriClip PRO2

Device	Size
PRO235	35 mm
PRO240	40 mm
PRO245	45 mm
PRO250	50 mm

**Selection Guide Device
CGG**

The AtriClip LAA Exclusion System is indicated for use in patients at high risk of thromboembolism for whom left atrial appendage exclusion is warranted.

Please review the Instructions for Use for a complete listing of contraindications, warnings, precautions and potential adverse events prior to using these devices.

1. Mean closure rate among the 3 studies is 98.1% per data on file at AtriCure, Inc. Ailawadi G. et al. J Thorac Cardiovasc Surg 2011;142:1002-9, 1009.e1; Gerdisch MW. et al. Innovations (Phila). 2022;17(6):463-470; Cartledge R et al., Interactive CardioVascular and Thoracic Surgery 34 (2022) 548-555

ATRICURE EUROPE B.V.

De Entree 260
1101 EE Amsterdam
Netherlands
Tel: +31 (0) 20-7005560
Fax: +31 (0) 20-7005561
SalesSupportEU@AtriCure.com
Europe.AtriCure.com

© AtriCure, Inc.
PM-EU-3094C-0928-G

AtriCure