


Hybrid AF™ Therapy EPI-Sense ST™ Coagulation Device

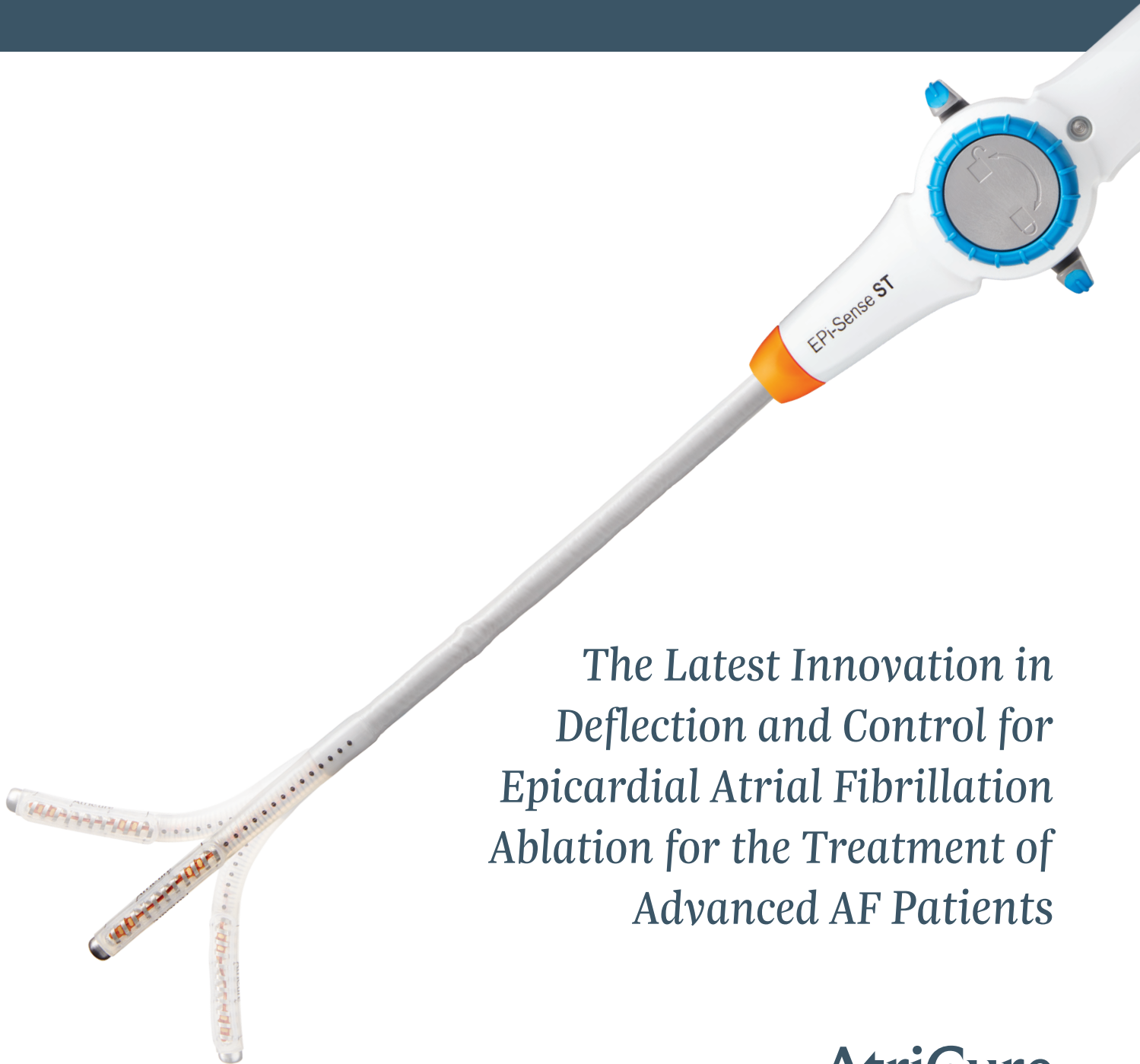

*The Latest Innovation in
Deflection and Control for
Epicardial Atrial Fibrillation
Ablation for the Treatment of
Advanced AF Patients*


AtriCure

Features

- Deflection of distal tip +/- 30° degrees
- Locking mechanism
- Torqueability (1:1 rotation)
- LED indicator light
- 3 cm device electrode length
- 2 distal and proximal sensing electrode pairs
- Irrigation/perfusion lumen
- Integrated suction

Benefits


- Steering lever allows the distal tip to deflect, allowing the operator to adjust the location of ablation coil
- Steering and braided catheter design provide control to position placement of the EPi-Sense ST
- Intuitive handle design and braided shaft allows for fine rotational and controlled positioning of the ablation coil
- LED on handle alerts user when device distal end exceeds temperature threshold

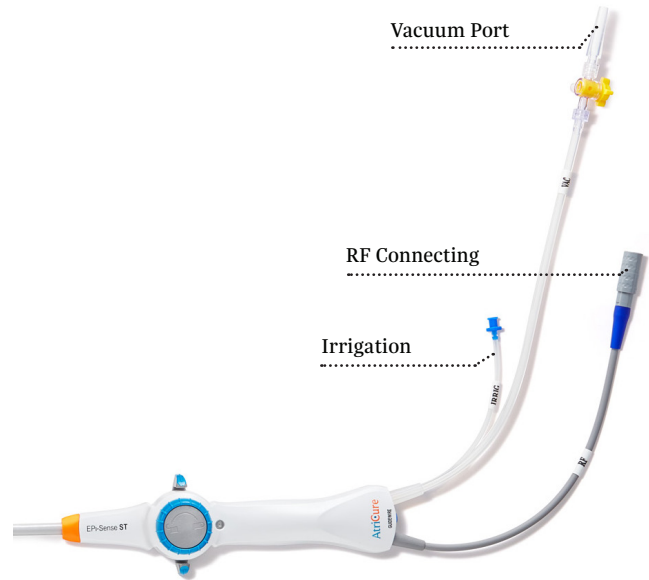

How EPi-Sense ST Works

Consistent Tissue Contact = Consistent Energy Transmission = Complete Lesions

Saline cooling solution cools the backside of the EPi-Sense ST device

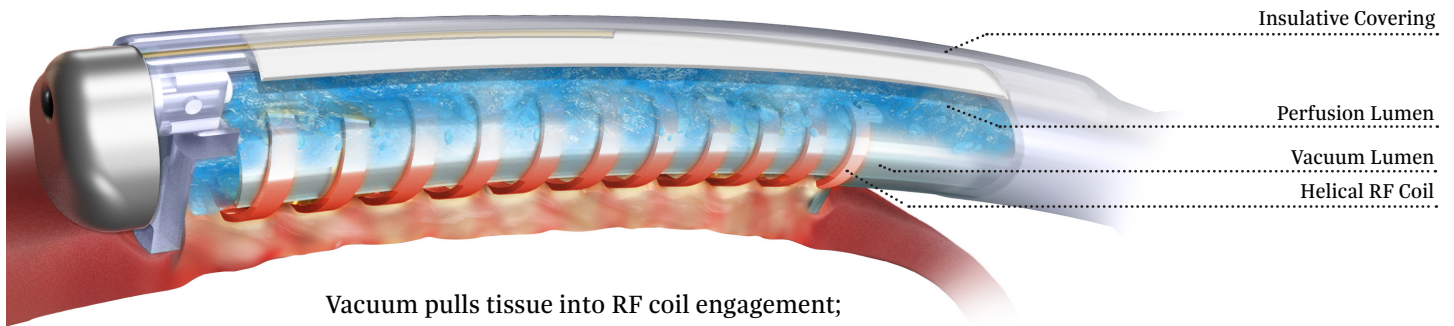

Vacuum pulls tissue into RF coil engagement;
perfusion over tissue conducts energy downward into
tissue while circulating blood absorbs excess heat.

EPi-Sense ST™ Coagulation Device

EU Indications: The EPi-Sense ST™ Coagulation Device is indicated for the epicardial treatment of arrhythmias including atrial fibrillation (AF) when augmented with an endocardial ablation, with the aim to restore normal sinus rhythm, reduce AF symptoms, and improve quality of life.

Please review the Instructions for Use for a complete listing of contraindications, warnings, precautions and potential adverse events prior to using this device.

ATRICURE EUROPE B.V.

De Entree 260
1101 EE
Amsterdam
The Netherlands
Tel. : +31 (0) 20-7005560
Fax : +31 (0) 20-7005561
SalesSupportEU@AtriCure.com
Europe.AtriCure.com

Epi-Sense ST Coagulation Device	
Device	Product Code
Epi-Sense ST Coagulation Device Kit w/30cm Cannula	EPIST-1
Epi-Sense ST Coagulation Device	EPIST
Epi-Sense ST RF Cable	CSK-2060
Cannula with Guide, 30cm	CSK-6131