




*cryo*SPHERE[®] MAX

60 seconds to manage
post-operative pain

AtriCure



cryoSPHERE MAX is designed to optimize the efficiency and efficacy of the Cryo Nerve Block procedure.

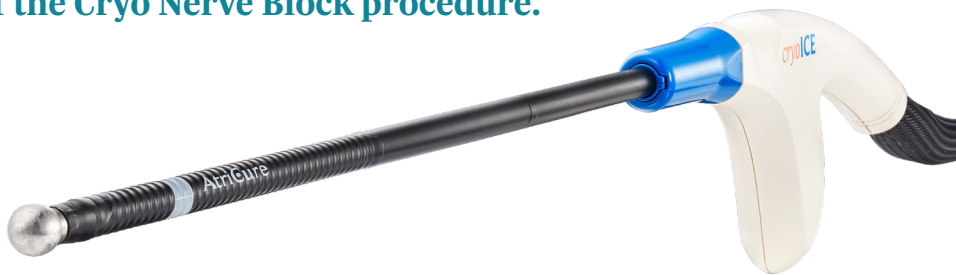

Increased Efficiency, Reduced Freeze Cycle

Featuring a 10 mm ball tip design, cryoSPHERE MAX provides improvements in heat extraction and ice ball formation.


cryoSPHERE MAX has 60% more tissue contact compared to cryoSPHERE’s 8 mm ball tip, optimizing the freeze cycle with a recommended freeze time of 60 seconds.

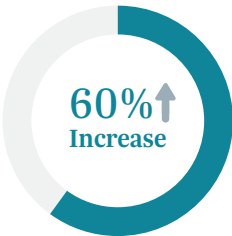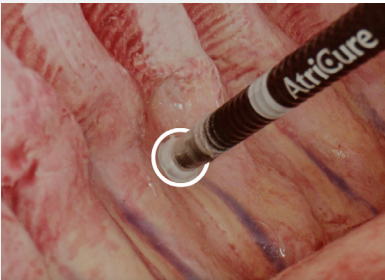

In surface area coverage⁺



In cryoablation freeze cycle⁺

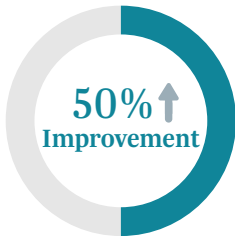


In nitrous oxide gas efficiency⁺

⁺ Compared to cryoSPHERE

[^] Reference on file at AtriCure

cryoSPHERE[®] MAX

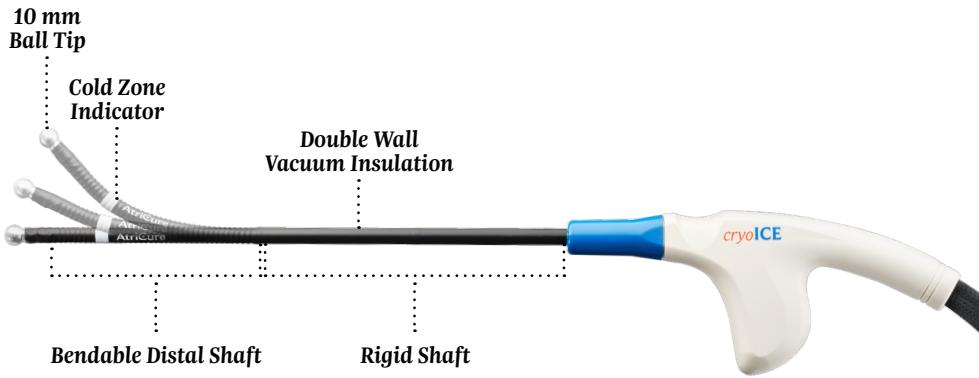

Features

Optimized Freeze Time	<ul style="list-style-type: none">• Recommended 60 second cryoablation cycle duration
10 mm Ball Tip	<ul style="list-style-type: none">• Larger surface area for maximizing heat extraction• Greater coverage for target tissues• 10 mm ball tip specifically engineered for cryoablation of peripheral nerves
Double Wall Vacuum Insulation	<ul style="list-style-type: none">• Designed to reduce risk of injury to non-target tissues that may contact the insulated shaft
Time To Therapeutic Temperature	<ul style="list-style-type: none">• Optimized time* to reach therapeutic temperature
Cold Zone Indicator	<ul style="list-style-type: none">• Provides visual reference that is approximately 4 cm* from the distal tip of the probe shaft**
Active Defrost	<ul style="list-style-type: none">• Allows the probe to be easily removed from the cryoablation site
Rigid Shaft	<ul style="list-style-type: none">• Stainless steel shaft provides rigidity and tactile feedback for user
Bendable Distal Shaft	<ul style="list-style-type: none">• Distal shaft bendable up to 140 degrees

*Reference on file at AtriCure

**Shaft distal of the Cold Zone Indicator can cause unintended cryoadhesion if touching unintended structures

Cryoablation Probes



11" equals 27,9 cm

17" equals 43,2 cm

cryoICE® BOX- V6

ACM2

cryoICE Cart & Upgrade Kits

Jansen Trolley (cart with double tank holder)

7.5-20 liter N₂O Cylinders, not provided by AtriCure

CMF1, Cryo Foot switch (optional)

EU Indications: The cryoICE® cryoSPHERE® MAX cryoablation probes are indicated for blocking pain associated with cardiac and/or thoracic procedures, by temporarily ablating intercostal peripheral nerves.

ATRICURE EUROPE B.V.

De Entrée 260

1101 EE Amsterdam Z.O.

Netherlands

Tel: +31 (0) 20-7005560

Fax: +31 (0) 20-7005561

SalesSupportEU@AtriCure.com

PM-EU-4925A-0227-G

AtriCure